

DTU





**DTU Health Technology
Bioinformatics**

Introduction to NGS technology

*Gabriel Renaud
Associate Professor
Section of Bioinformatics
Technical University of Denmark
gabriel.reno@gmail.com*

Outline

- 2nd generation NGS
- Illumina movietime!
- Your turn to basecall
- 3rd generation NGS

2 main types of approaches

1) Amplify and sequence one base at a time

1:A 2:G 3:G 4:T = AGGT

2) Amplify and count how many of the same base

1:1A 2:2G 3:1T = AGGT

Second generation sequencing

- Illumina sits on 75% of the market



Illumina




BGISEQ



454




Ion Torrent

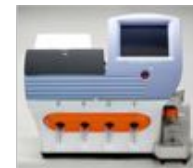


General library preparation steps

1. Create library molecules
2. Amplification (PCR)
3. Massive parallel sequencing (strength over Sanger)

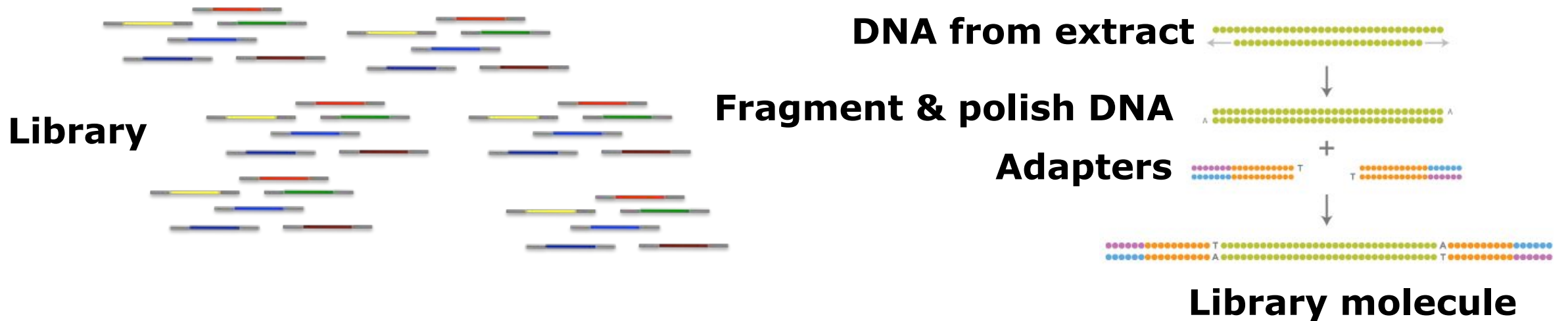


Amplification and immobilization

- Emulsion PCR (454, SOLiD, IonTorrent): Water, oil, beads, one DNA template/droplet

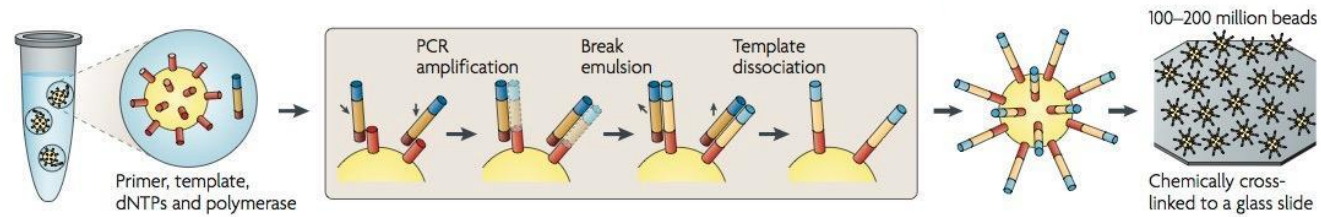


Amplification and immobilization

- Emulsion PCR (454, SOLiD, IonTorrent): Water, oil, beads, one DNA template/droplet

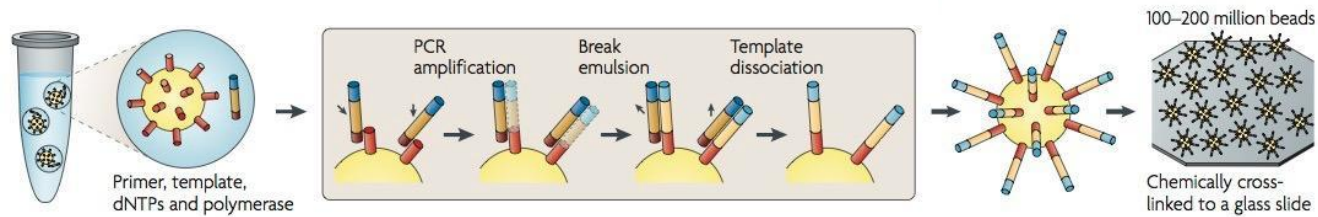

Bridge PCR (Illumina): One DNA template/cluster, primers on surface, grow by bridging primers

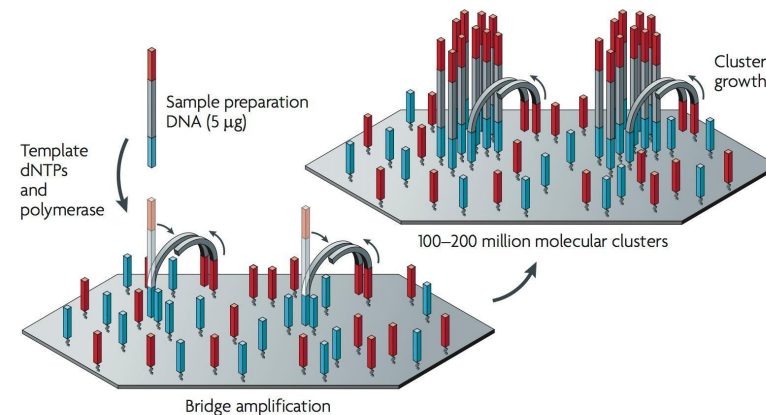

2 main types of approaches

1) Amplify and sequence one base at a time

1:A 2:G 3:G 4:T = AGGT

2) Amplify and count how many of the same base

1:1A 2:2G 3:1T = AGGT


Illumina sequencing

corporate propaganda:

<https://www.youtube.com/watch?v=HMyCqWhwB8E>

Amplicon sequencing on Illumina

- Why can't you just fill your Illumina flow cell with amplicon libraries (i.e. the same sequence over and over)?

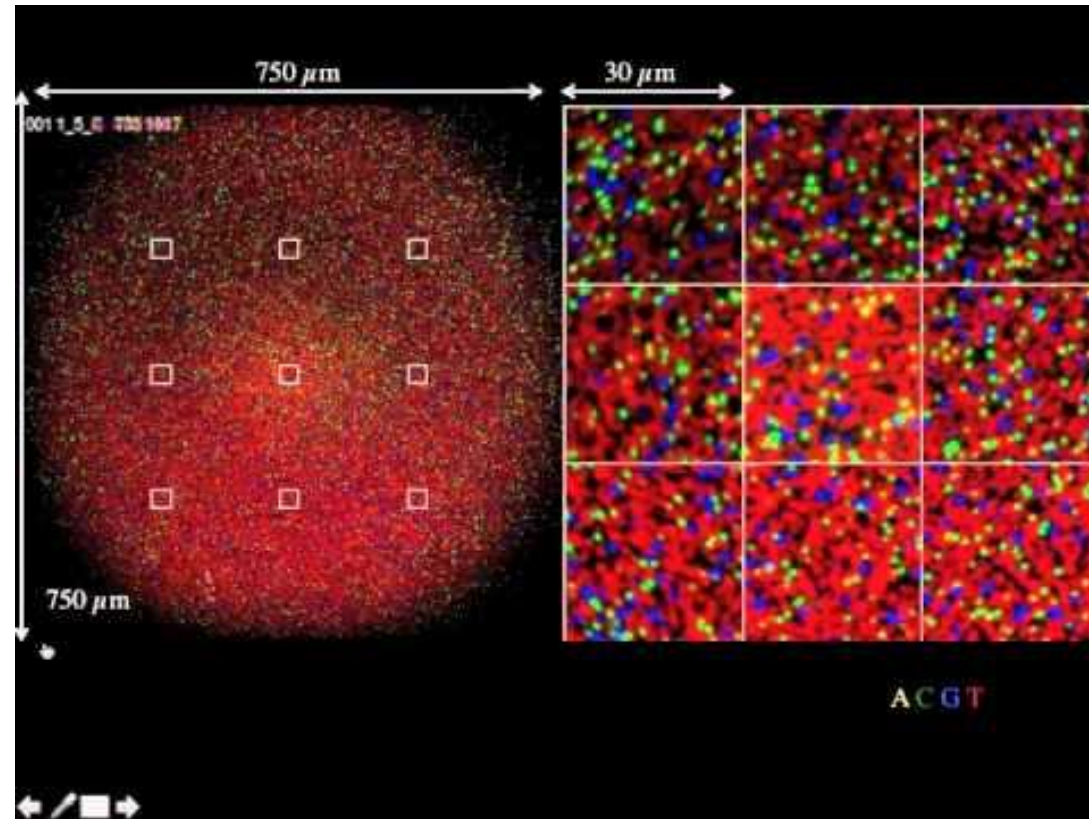


Fluorescence detection

Illumina - Cyclic reversible termination

Add all dNTPs labelled w.
diff dye

Create four-color image

Cleave dye and repeat next
cycle

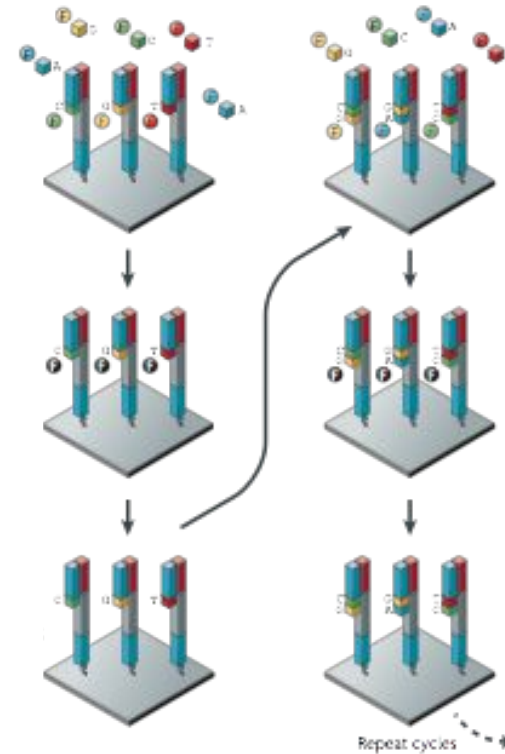

2G: Imaging handout



Illumina 1: _____

Illumina 2: _____

—
—
—

2G: Imaging handout Answers!



Illumina 1: _____

Illumina 2: _____


TOP: **CATCGT**


BOTTOM: **CCCCCC**

—
—
—

Illumina: Quality deterioration


- Quality goes down
- Especially 2nd read
- Can you think of why?

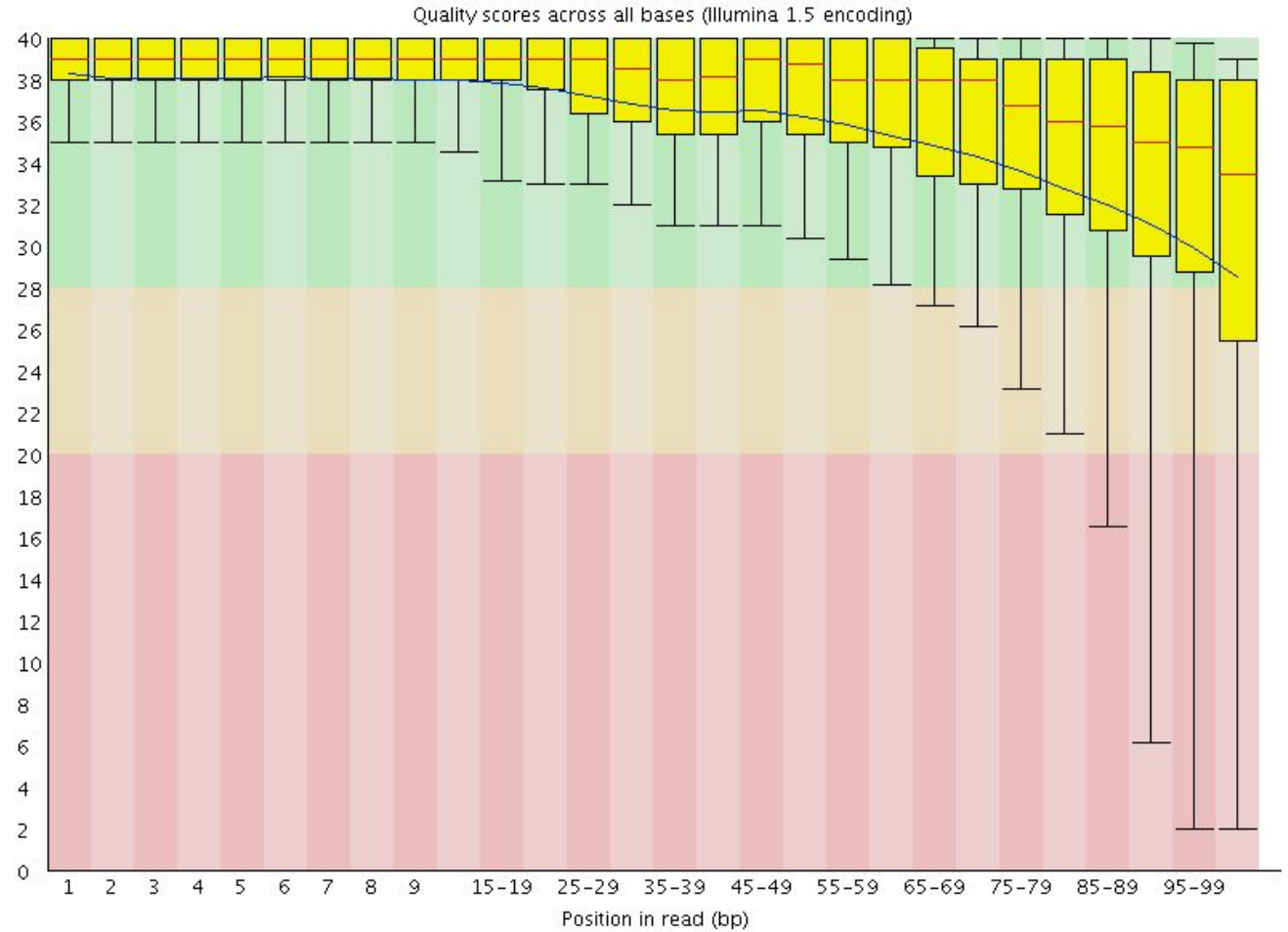

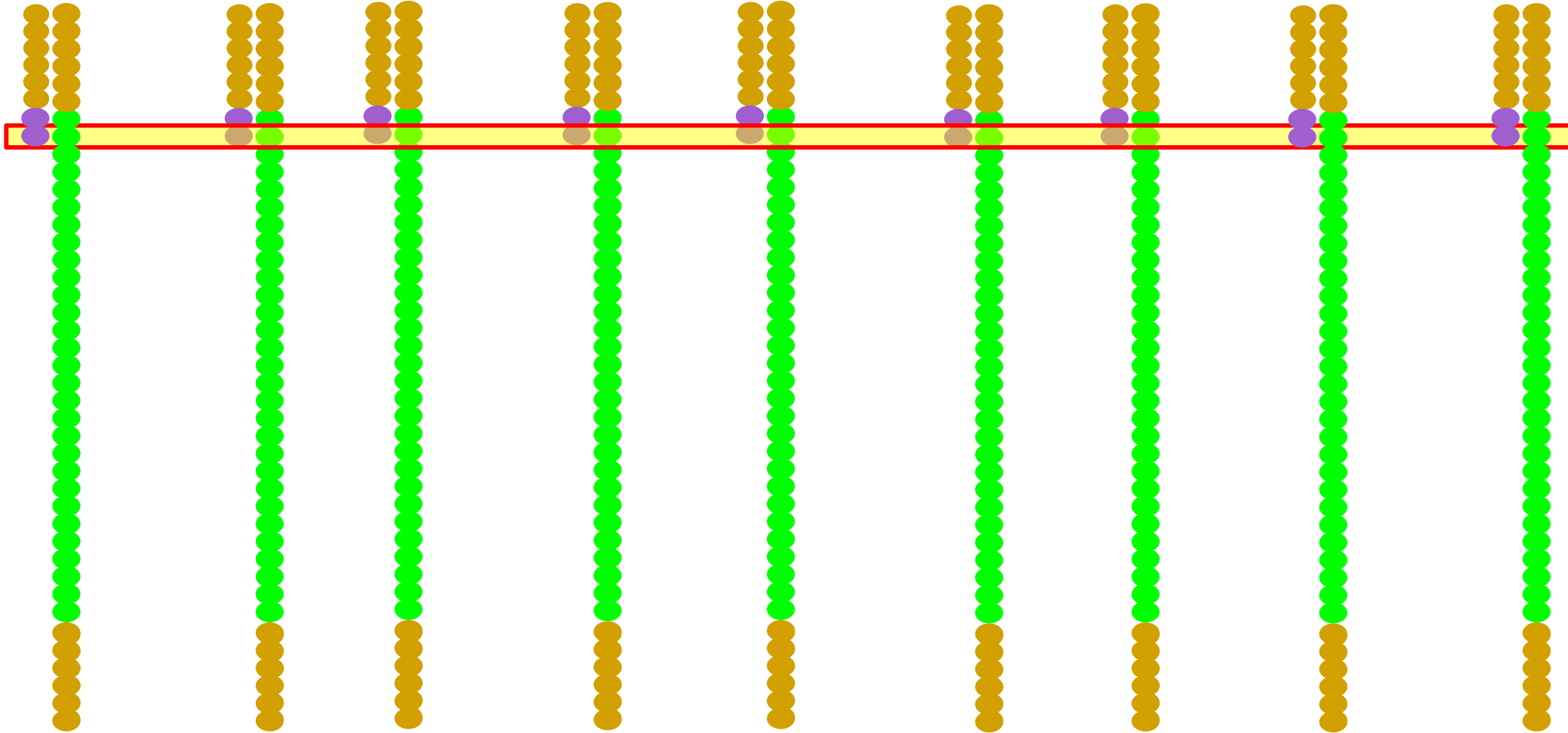

Phasing


Prephasing

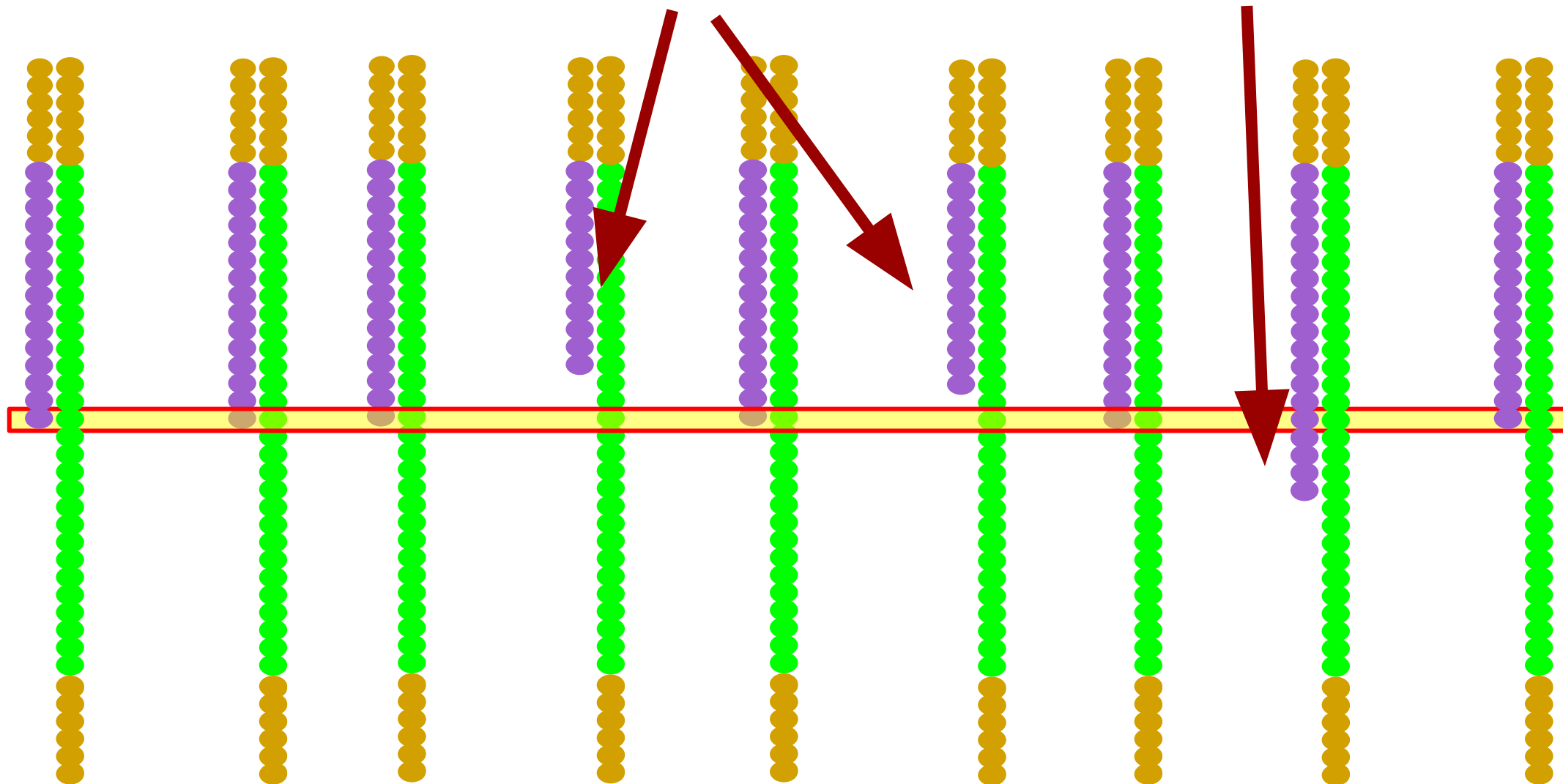

Illumina: Quality deterioration

- Quality goes down
- Especially 2nd read
- Can you think of why?

- Efficiency of incorporation
- Phasing
- Prephasing

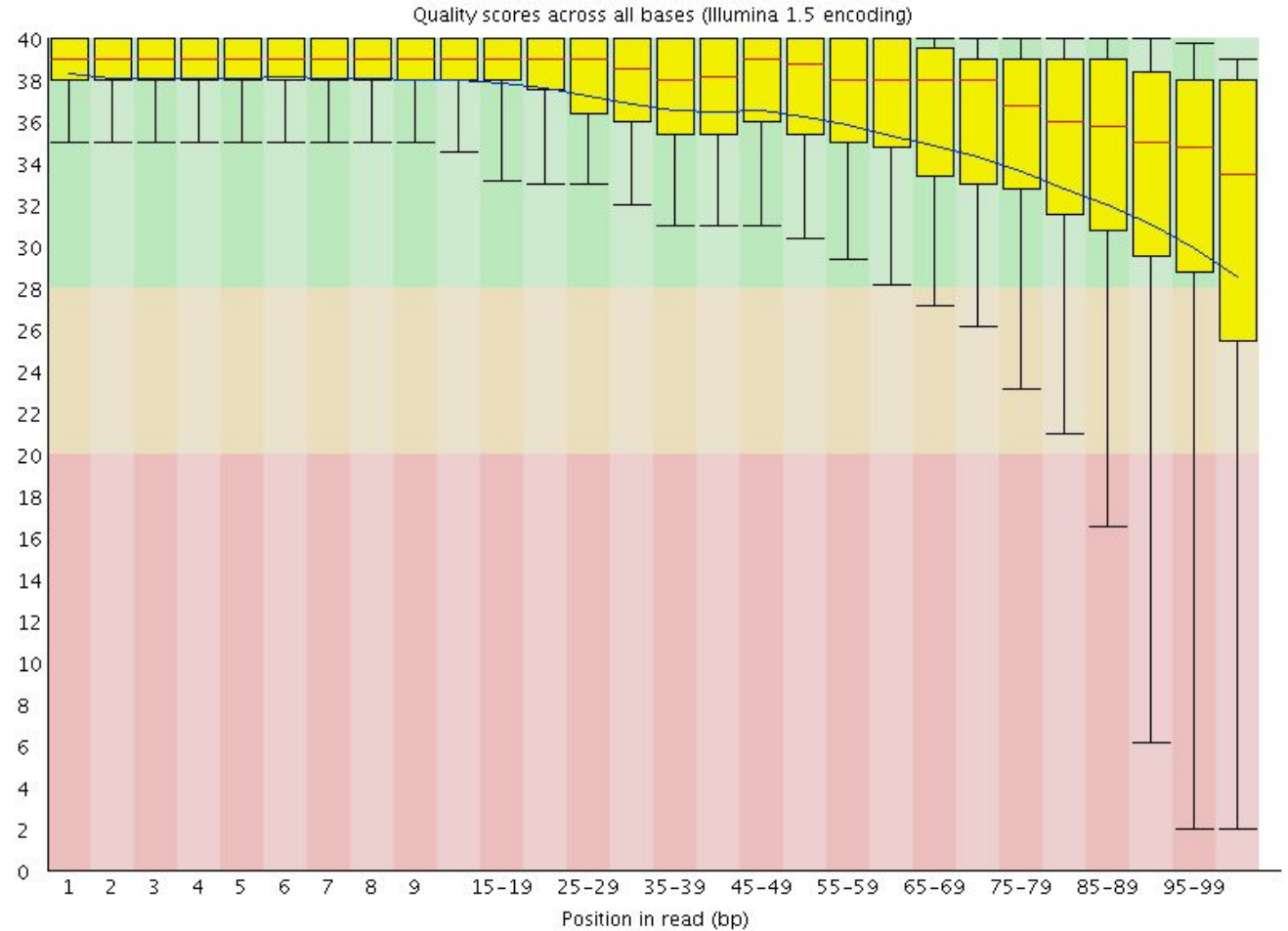

Brief side note about multiplexing/demultiplexing


- If we sequence a small virus (ex: bacteriophage Phi-X174 with a genome size of 5386 nucleotides), do we need 1B reads?
- Idea to save costs: pool multiple samples together on the same run

Brief side note about spike-in

- How to know if the sequencing run was successful (low error rate)?
- Idea: Let's spike-in a small virus (ex: bacteriophage Phi-X174 with a genome size of 5386 nucleotides)


NextSeq/NovaSeq (2015-)

- Chemistry is not based 4 dyes (as before) but 2 dyes
 - T (red), C (green), A (both) and G (none = “dark”)
 - Faster processing rate and cheaper reagents
 - Slightly increases error rate
 - Problem with G stretches because G is not dyed

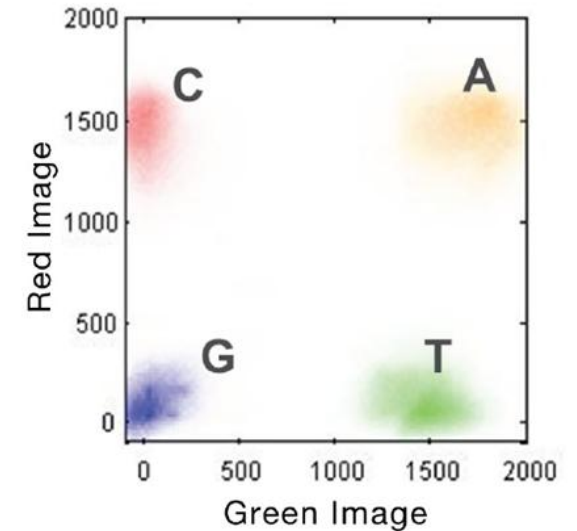

source: Illumina

1 dye, 2 images

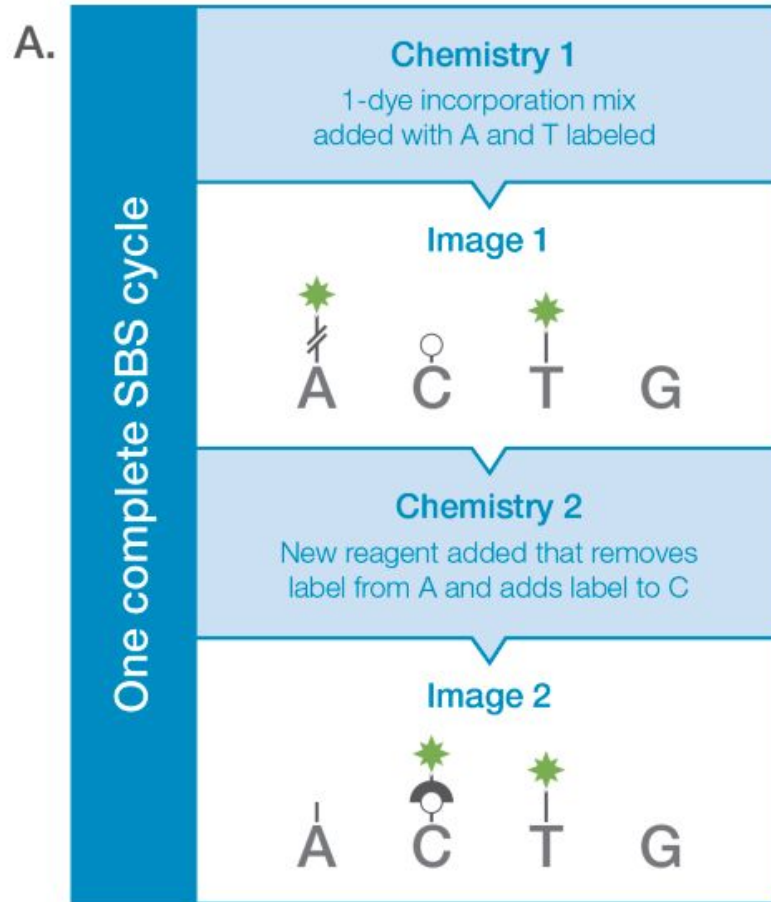


B.

Image 1	Image 2	Result
ON	OFF	A
OFF	ON	C
ON	ON	T
OFF	OFF	G

source: Illumina

Patterned flowcell


- Patterned wells
- Solves overloading flowcell
- More duplicates

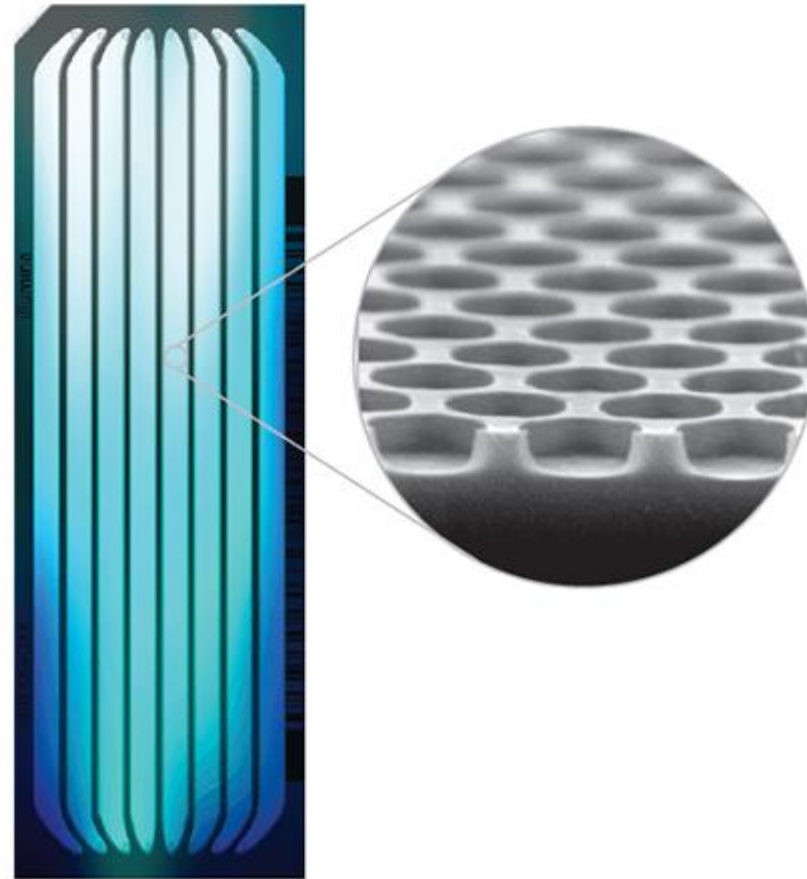


source: Illumina

BGI-Seq


ssDNA -> DNA nanoballs

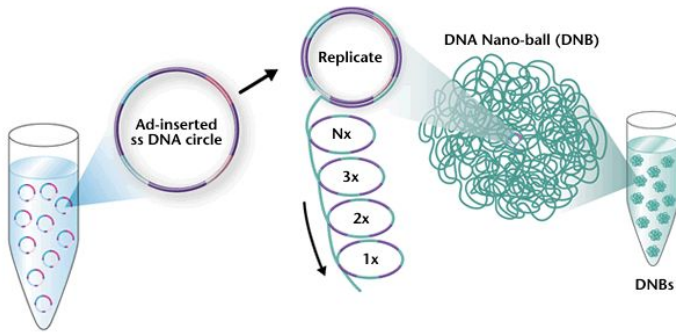


Place DNBs into each spot

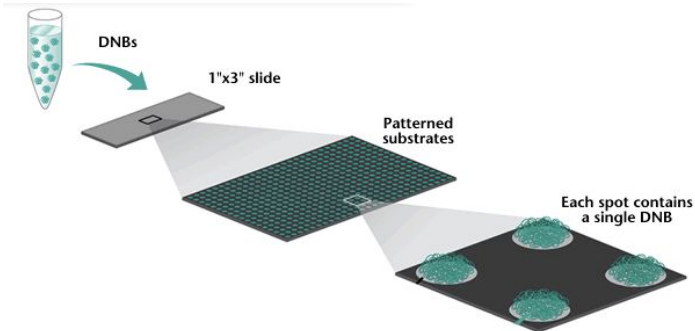

Use silicon chips with sticky spots

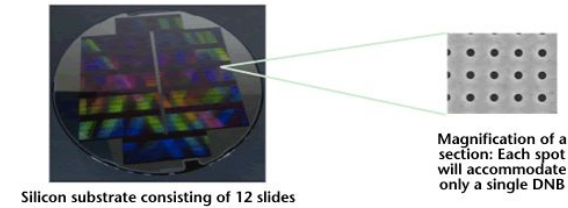

Sequence using ligase and fluorescent labeled probes

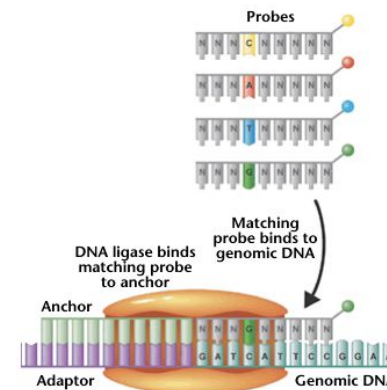

2020

PLOS ONE

RESEARCH ARTICLE

Comparative analysis of novel MGISEQ-2000 sequencing platform vs Illumina HiSeq 2500 for whole-genome sequencing

Dmitriy Korostin¹, Nikolay Kulemin^{1,2}, Vladimir Naumov², Vera Belova^{1*}, Dmitriy Kwon³, Alexey Gorbachev³

¹ Pirogov Russian National Research Medical University, Moscow, Russia, ² Zenome.io, Ltd., Moscow, Russia, ³ Company Helicon, Ltd., Moscow, Russia

* verusik.belova@gmail.com

Abstract

The MGISEQ-2000 developed by MGI Tech Co. Ltd. (a subsidiary of the BGI Group) is a new competitor of such next-generation sequencing platforms as NovaSeq and HiSeq (Illumina). Its sequencing principle is based on the DNB and the cPAS technologies, which were also used in the previous version of the BGISeq-500 device. However, the reagents for MGISEQ-2000 have been refined and the platform utilizes updated software. The cPAS method is an advanced technology based on the cPAL previously created by Complete Genomics. In this paper, the authors compare the results of the whole-genome sequencing of a DNA sample from a Russian female donor performed on MGISEQ-2000 and Illumina HiSeq 2500 (both PE150). Two platforms were compared in terms of sequencing quality, number of errors and performance. Additionally, we performed variant calling using four different software packages: Samtools mpileup, Strelka2, Sentieon, and GATK. The accuracy of SNP detection was similar in the data generated by MGISEQ-2000 and HiSeq 2500, which was used as a reference. At the same time, a separate indel analysis of the overall error rate revealed similar FPR values and lower sensitivity. It may be concluded with confidence that the data generated by the analyzed sequencing systems is characterized by comparable magnitudes of error and that MGISEQ-2000 and HiSeq 2500 can be used interchangeably for similar tasks like whole genome sequencing.

Comparative Performance of the MGISEQ-2000 and Illumina X-Ten Sequencing Platforms for Paleogenomics

Kongyang Zhu^{1†}, Panxin Du^{2†}, Jianxue Xiong², Xiaoying Ren³, Chang Sun², Yichen Tao², Yi Ding³, Yiran Xu², Hailiang Meng², Chuan-Chao Wang^{1*} and Shao-Qing Wen^{2,3*}

¹State Key Laboratory of Cellular Stress Biology, School of Life Sciences, State Key Laboratory of Marine Environmental Science, Department of Anthropology and Ethnology, Institute of Anthropology, School of Sociology and Anthropology, Xiamen University, Xiamen, China, ²MCE Key Laboratory of Contemporary Anthropology, Department of Anthropology and Human Genetics, School of Life Sciences, Fudan University, Shanghai, China, ³Institute of Archaeological Science, Fudan University, Shanghai, China

The MGISEQ-2000 sequencer is widely used in various omics studies, but the performance of this platform for paleogenomics has not been evaluated. We here compare the performance of MGISEQ-2000 with the Illumina X-Ten on ancient human DNA using four samples from 1750 BCE to 60 CE. We found there were only slight differences between the two platforms in most parameters (duplication rate, sequencing bias, θ , δS , and λ). MGISEQ-2000 performed well on endogenous rate and library complexity although X-Ten had a higher average base quality and lower error rate. Our results suggest that MGISEQ-2000 and X-Ten have comparable performance, and MGISEQ-2000 can be an alternative platform for paleogenomics sequencing.

2021



DATA NOTE

Comparative analysis of 7 short-read sequencing platforms using the Korean Reference Genome: MGI and Illumina sequencing benchmark for whole-genome sequencing

Hak-Min Kim¹, Sungwon Jeon^{2,3}, Oksung Chung¹, Je Hoon Jun¹, Hui-Su Kim², Asta Blazyte^{2,3}, Hwang-Yeol Lee¹, Youngseok Yu¹, Yun Sung Cho¹, Dan M. Bolser^{4,*} and Jong Bhak^{1,2,3,4,5,*}

¹Clinomics Inc., Ulsan National Institute of Science and Technology (UNIST), UNIST-gil 50, Eonyang-eup, Ulju-gun, Ulsan, 44919, Republic of Korea; ²Korean Genomics Center (KOGIC), Ulsan National Institute of Science and Technology (UNIST), UNIST-gil 50, Eonyang-eup, Ulju-gun, Ulsan, 44919, Republic of Korea; ³Department of Biomedical Engineering, School of Life Sciences, Ulsan National Institute of Science and Technology (UNIST), UNIST-gil 50, Eonyang-eup, Ulju-gun, Ulsan, 44919, Republic of Korea; ⁴Geromics Ltd., 222 Mill Road, Cambridge, CB1 3NF, United Kingdom and ⁵Personal Genomics Institute (PGI), Genome Research Foundation, Osong saengmyong1ro, Cheongju, 28160, Republic of Korea


GigaScience, 10, 2021, 1-9

doi: 10.1093/gigascience/gjab014
Data Note

Conclusion: BGI = Illumina in terms of errors but cheaper

BGI-Seq

Google scholar articles on BGISeq per year

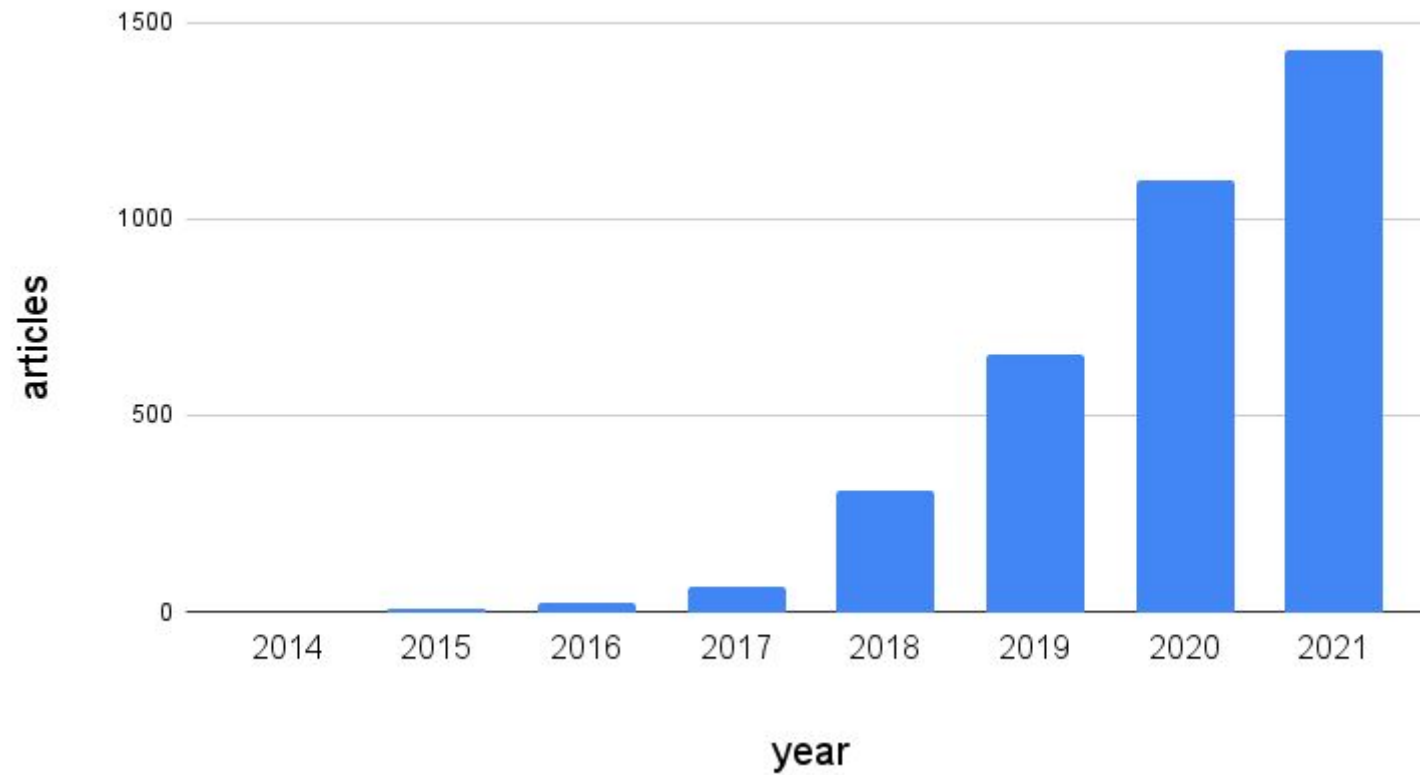

2 main types of approaches

1) Amplify and sequence one base at a time


1:A 2:G 3:G 4:T = AGGT

2) Amplify and count how many of the same base

1:1A 2:2G 3:1T = AGGT

454: Pyrosequencing

1. Load template beads into wells
2. Flow one dNTP across wells
3. Polymerase incorporates nucleotide
4. Release of PP_i leads to light
5. Light intensity= # of bases
6. Imaging, next dNTP

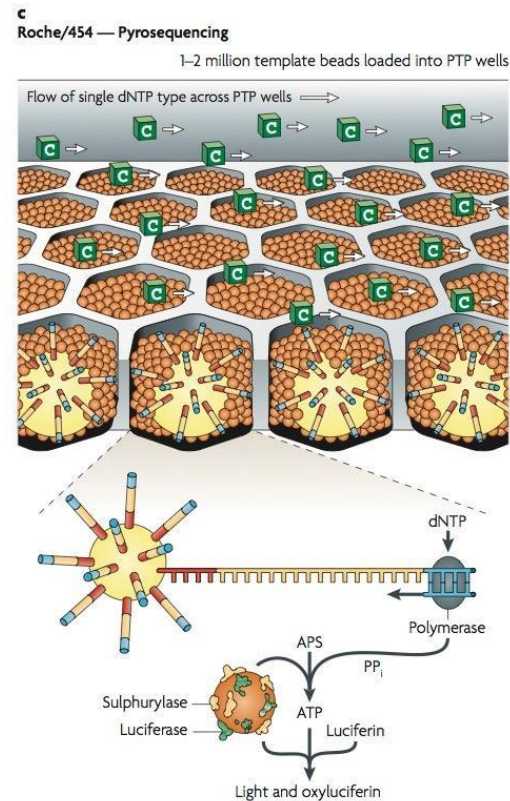

2G: Imaging handout




454: _____

-

-

2G: Imaging handout Answers!


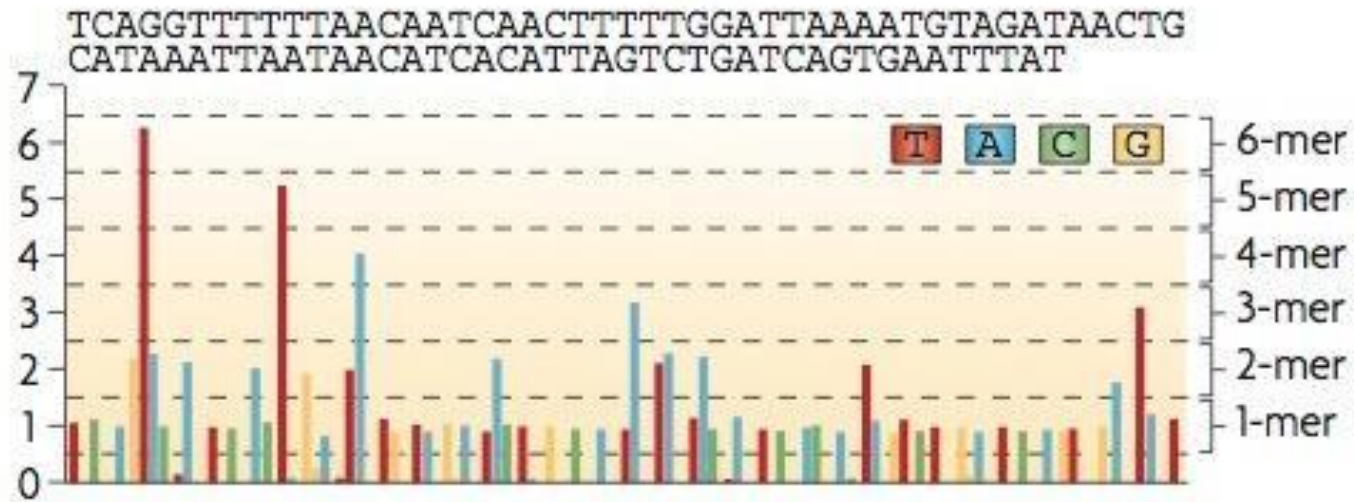


454: _____

-
 -

Ion Torrent

- Corporate propaganda: <https://www.youtube.com/watch?v=zBPKj0mMcDg>
- Similar principle to 454
- Library: Emulsion PCR
- Based on semiconductors
- Detection is based on H ions (pH) changes

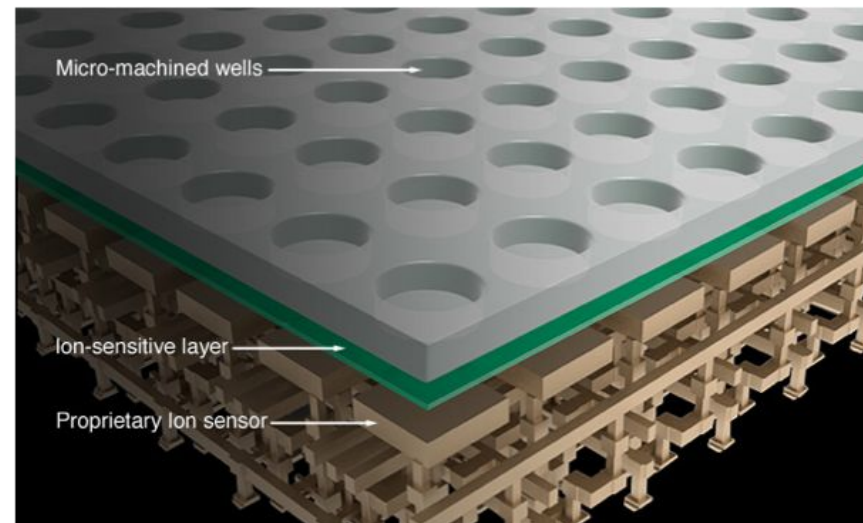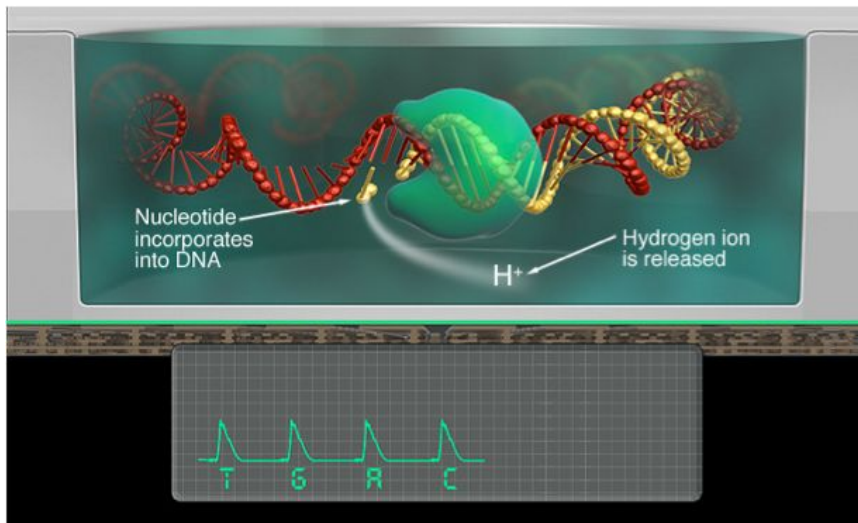

Let's remember the types of errors

mismatch

AGCAATCTCAATTACAAATATACACCAACAAA

AGCAATCTCAATTACAGATATACACCAACAAA

insert

AGCAATCTCAATTACA-AATATACACCAACAA

AGCAATCTCAATTACACAATATACACCAACAA

deletion

AGCAATCTCAATTACAAATATACACCAACAA

AGCAATCTCAATTACA-ATATACACCAACAA

Quiz!

Which of the the 2 main types of approaches would be more prone to indels?

1) Amplify and sequence one base at a time

1:A 2:G 3:G 4:T = AGGT


2) Amplify and count how many of the same base

1:1A 2:2G 3:1T = AGGT

Technology	read length	# of reads	errors?
Sanger	400 to 900 bp	96	mm 0.01%
Illumina MiSeq	2x 200-300bp	20-30 M per flow cell	mm 0.1-0.2%
Illumina HiSeq	2x 100-150bp	~2G per flow cell	mm 0.1-0.2%
Illumina NovaSeq	2x 100-250bp	~6-10G per flow cell	mm 0.1%?
DNBSEQ-T7	2x 100-200bp	~5-10G per flow cell	<mm 0.1%
Ion Torrent	~200-400 bp	5-50M reads	indel 0.46 to 2.4%

3rd generation

1977 1980 1983 1986 1989 1992 1995 1998 2001 2004 2007 2010 2013 2016 2019

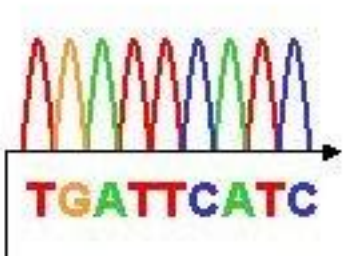

Sanger



454



PacBio



Ion Torrent



Oxford Nanopore

SOLiD



BGI



Illumina

3rd generation

- Single-molecule sequencing
- No amplification -> less bias -> observations are more independent



Helicos



PacBio




Oxford Nanopore

Oxford Nanopore


- Literal nanopores
- Current per base
- Non-random errors
- <https://www.youtube.com/watch?v=RcP85JHLmnl>
- Very high error rate

“If a nanopore was the size of a fist, a 1MB strand of DNA passing through that nanopore would be 2 miles (3.2 km) long”
 -Adam Philippy, NHGRI

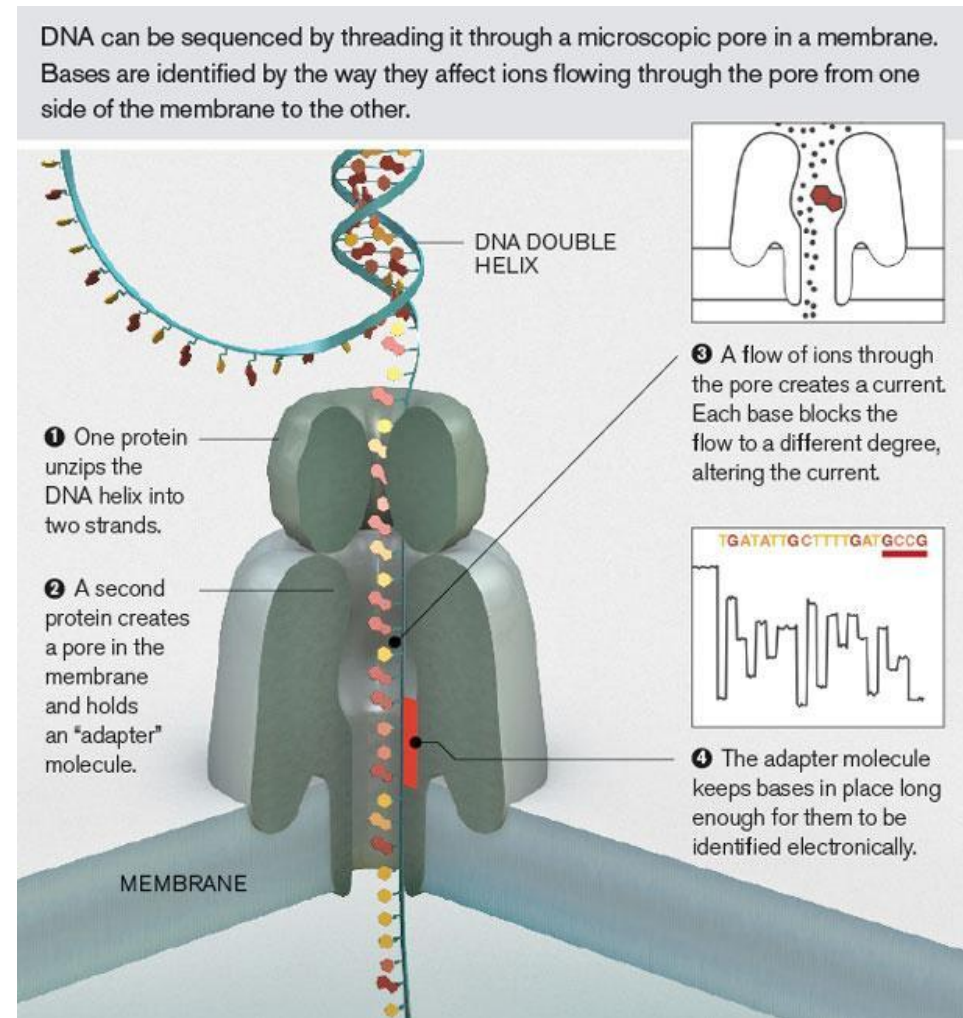

Oxford Nanopore

- Hairpin allows double sequencing (2D)

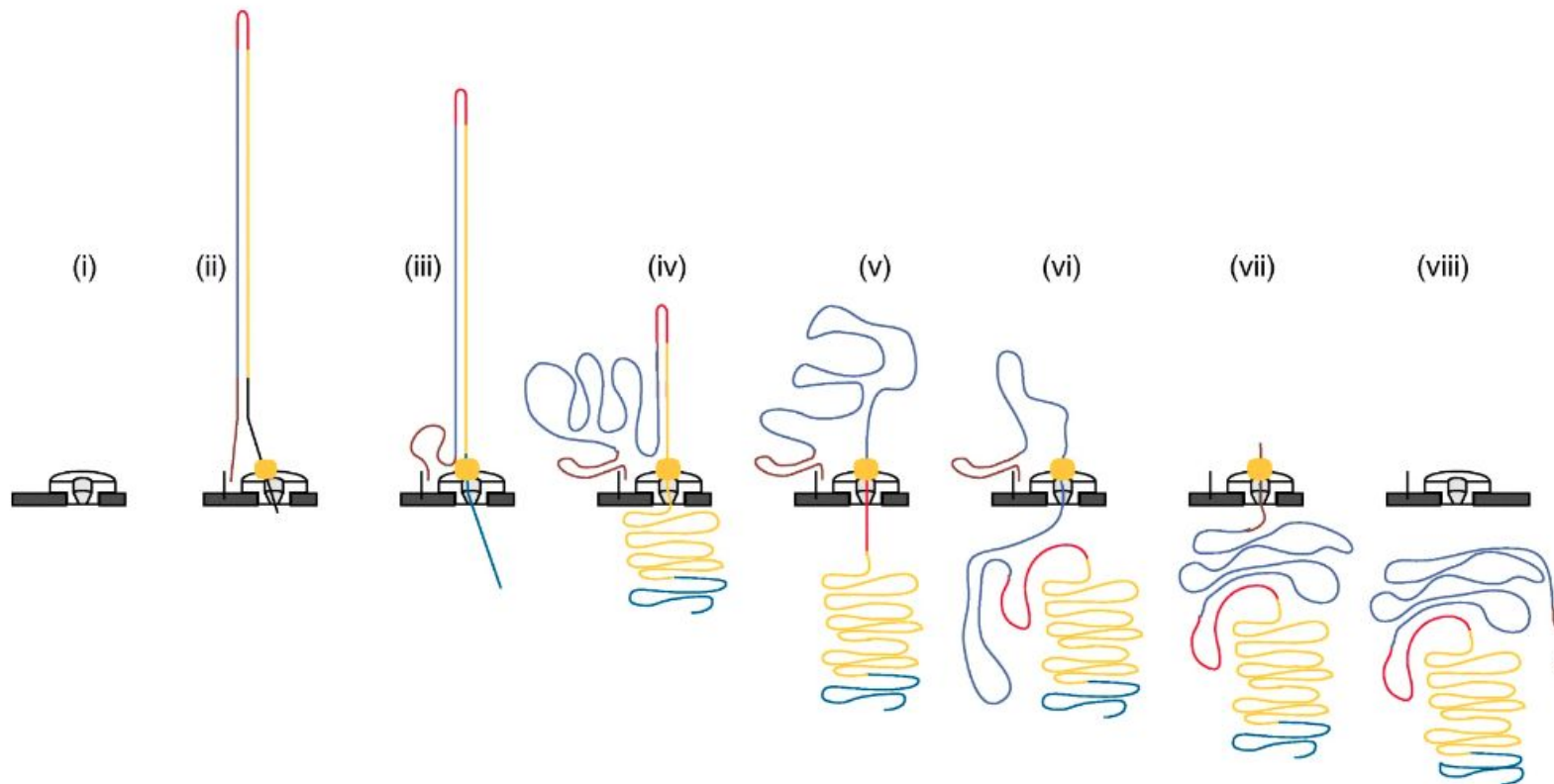

Jain, M., Olsen, H.E., Paten, B. et al. The Oxford Nanopore MinION: delivery of nanopore sequencing to the genomics community. *Genome Biol* 17, 239 (2016). <https://doi.org/10.1186/s13059-016-1103-0>

Cheap & mobile

- Long reads, low quality
- Low establishment and maintenance costs
- Portability



PacBio: Single-molecule real-time (SMRT) sequencing

- Expensive machinery
- Not very portable



PacBio


- Flexibility
 - Long but low quality or shorter but better reads
 - Robust
 - https://www.youtube.com/watch?v=_ID8JyAbwEo
 - New 2019: HiFi read same fragment multiple times

High-throughput sequencing



Library preparation

SMRTBell 'template'

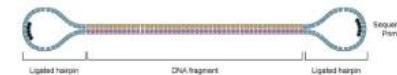

Standard 'Sequencing'



Circular 'Consensus' Sequencing'



Single pass



Continued generations of reads

Multiple passes




Tiny wells

- 1 million wells per cell
- Hit the lights

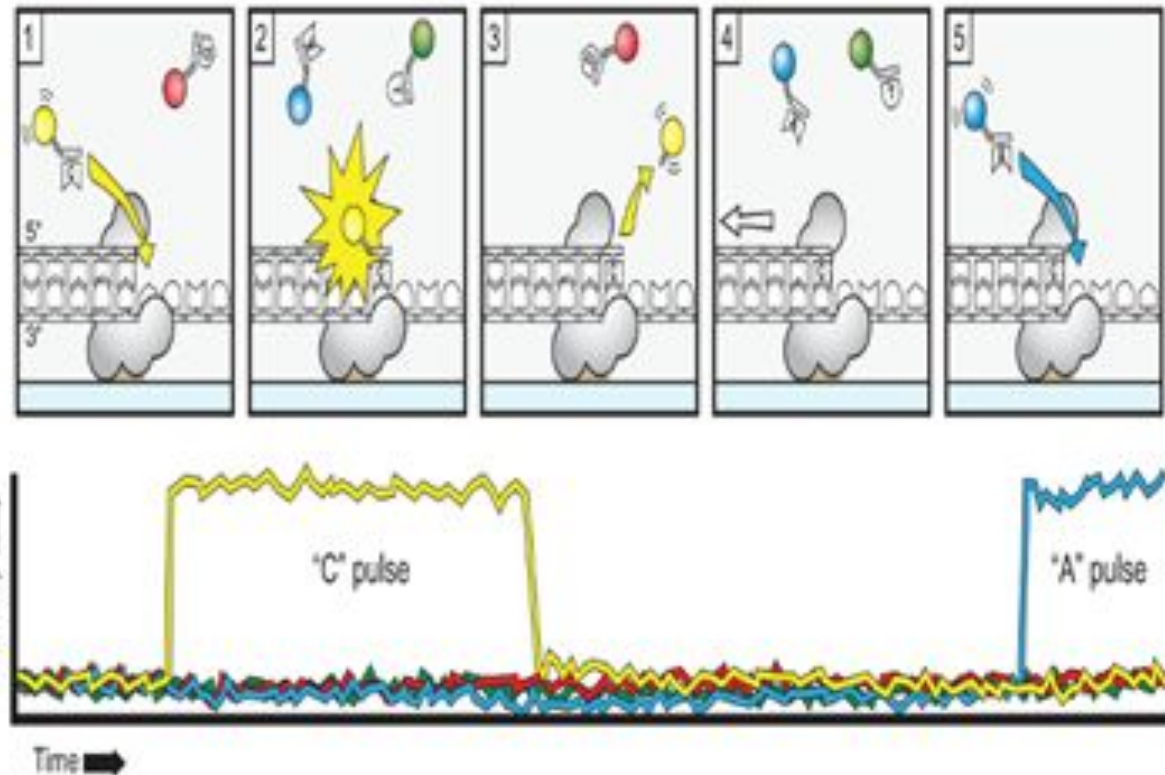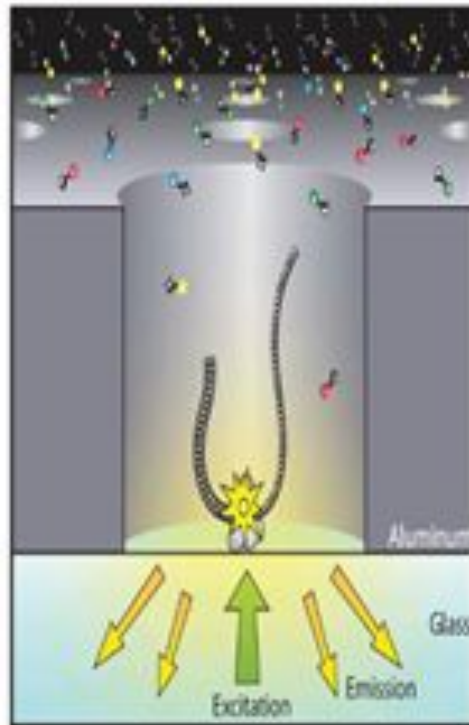


Technology	read length	# of reads	errors?
Oxford Nanopore	avg. 2 kbp-10 kbp	2M-6G	2021 update: ~3% 1D: indel+mm 20% 2D: indel+mm 7%
PacBio	10-20 kbp	500k-4M	indel+mm 13-15% HiFi: indel 1%+mm 0.1%

The Oxford Nanopore error rates are already outdated...



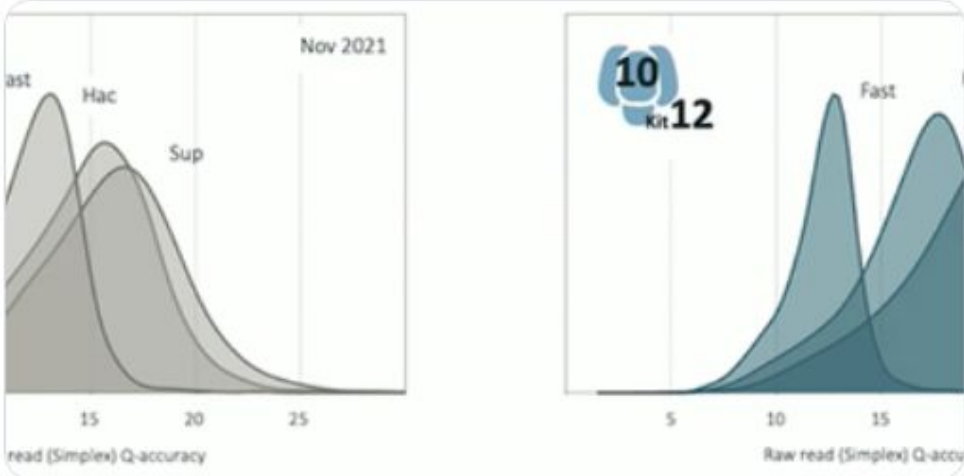

 **Oxford Nanopore** ✓
@nanopore

Kit 12 — which includes the latest “Q20+” chemistry for >Q20 (>99%) raw read and around Q30 (99.9%) Duplex sequencing accuracy — and flow cells containing the latest pore (R10.4), are now available to buy in Store: [store.nanoporetech.com/sample-prep.ht...](https://store.nanoporetech.com/sample-prep.html) #nanoporeconf



9:27 PM · Dec 1, 2021 · Twitter Web App

17 Retweets **2** Quote Tweets **34** Likes

Article | [Published: 09 September 2021](#)


Performance assessment of DNA sequencing platforms in the ABRF Next-Generation Sequencing Study

[Jonathan Foox](#), [Scott W. Tighe](#), [...] [Christopher E. Mason](#) 

[Nature Biotechnology](#) **39**, 1129–1140 (2021) | [Cite this article](#)

5529 Accesses | 171 Altmetric | [Metrics](#)

 An [Author Correction](#) to this article was published on 11 October 2021

 This article has been [updated](#)

Abstract

Assessing the reproducibility, accuracy and utility of massively parallel DNA sequencing


Takeaways:

Short reads

- Illumina cheapest
- BGI most accurate

Long reads:

- Most mapping with PacBio
- Oxford/Pacbio good with repeats

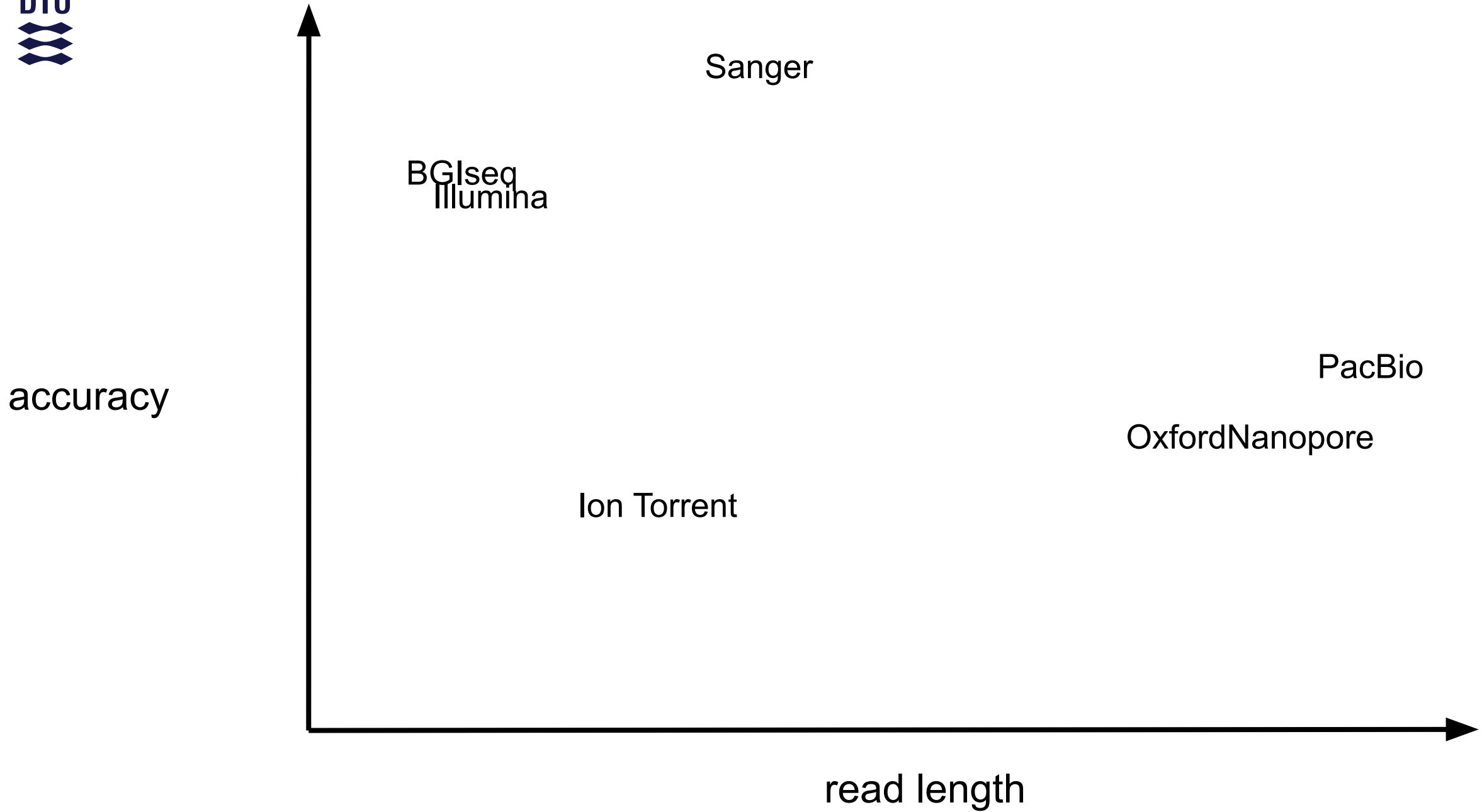

font size = throughput

accuracy

BGIseq
Illumina

that is Sanger

Ion Torrent

OxfordNanopore

PacBio

read length

Summary

- I did not mention a very important factor: **cost**
- Each tech has advantages, pick the most appropriate for your question
- Illumina is the current workhorse
 - Great for many applications
- Long read technology
 - Adding information
 - Resolves difficult regions during genome assembly

Article | [Open Access](#) | Published: 14 July 2020

Telomere-to-telomere assembly of a complete human X chromosome

Karen H. Miga , Sergey Koren, Arang Rhie, Mitchell R. Vollger, Ariel Gershman, Andrey Bzikadze, Shelise Brooks, Edmund Howe, David Porubsky, Glennis A. Logsdon, Valerie A. Schneider, Tamara Potapova, Jonathan Wood, William Chow, Joel Armstrong, Jeanne Fredrickson, Evgenia Pak, Kristof Tigyi, Milinn Kremitzki, Christopher Markovic, Valerie Maduro, Amalia Dutra, Gerard G. Bouffard, Alexander M. Chang, Nancy F. Hansen, Amy B. Wilfert, Françoise Thibaud-Nissen, Anthony D. Schmitt, Jon-Matthew Belton, Siddarth Selvaraj, Megan Y. Dennis, Daniela C. Soto, Ruta Sahasrabudhe, Gulhan Kaya, Josh Quick, Nicholas J. Loman, Nadine Holmes, Matthew Loose, Urvashi Surti, Rosa ana Risques, Tina A. Graves, Lindsay, Robert Fulton, Ira Hall, Benedict Paten, Kerstin Howe, Winston Timp, Alice Young, James C.